




Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

● VT5 Set ○ VT8 Set ● VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs																		
Implant Direct®			Legacy™ I															
			Soft Bone							Hard Bone (Mandible)								
										In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. Click here to view PDF.								
Geometry	Major Ø	Minor Ø	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Densah® Bur Block Display	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Densah® Bur Block Display	
Taper	3.7	3.1	Pilot	VT1525 (2.0)	VT2535* (3.0)	—	—	—		Pilot	VT1525 (2.0)	VT1828 (2.3)	VT2535* (3.0)	—	—	—	—	
Taper	4.7	4.1	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545* (4.0)	—	—		Pilot	VT1525 (2.0)	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545* (4.0)	—	—	
Taper	5.7	5.1	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545 (4.0)	VT4555* (5.0)	—		Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4555* (5.0)	—	

*Denotes implant placement. *Clinician experience and judgement should be used in conjunction with the Densifying Reference Guide recommendation and suggested use protocols. W042 REV04

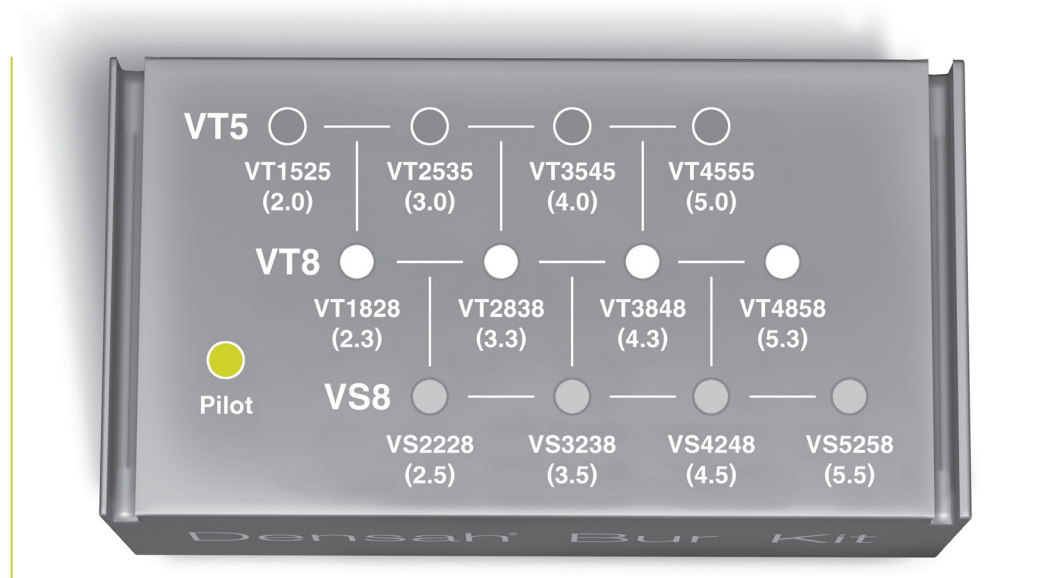


Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

- VT5 Set
- VT8 Set
- VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs																	
Implant Direct®			Legacy™ 2, 3, 4														
			Soft Bone							Hard Bone (Mandible)							
										In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. Click here to view PDF.							
Geometry	Major Ø	Minor Ø	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Densah® Bur Block Display	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Densah® Bur Block Display
Taper	3.2	2.6	Pilot	VT1828 (2.3)	VS2228* (2.5)	—	—	—		Pilot	VT1828 (2.3)	VS2228* (2.5)	—	—	—	—	
Taper	3.7	3.1	Pilot	VT1525 (2.0)	VT2535* (3.0)	—	—	—		Pilot	VT1525 (2.0)	VT1828 (2.3)	VT2535* (3.0)	—	—	—	
Taper	4.2	3.6	Pilot	VT1828 (2.3)	VT2838* (3.3)	—	—	—		Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VS3238* (3.5)	—	—	
Taper	4.7	4.1	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545* (4.0)	—	—		Pilot	VT1525 (2.0)	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545* (4.0)	—	
Taper	5.2	4.6	Pilot	VT1828 (2.3)	VT2838 (3.3)	VT3848* (4.3)	—	—		Pilot	VT1828 (2.3)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VS4248* (4.5)	
Taper	5.7	5.1	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545 (4.0)	VT4555* (5.0)	—		Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4555* (5.0)	

*Denotes implant placement. *Clinician experience and judgement should be used in conjunction with the Densifying Reference Guide recommendation and suggested use protocols. W042 REV04

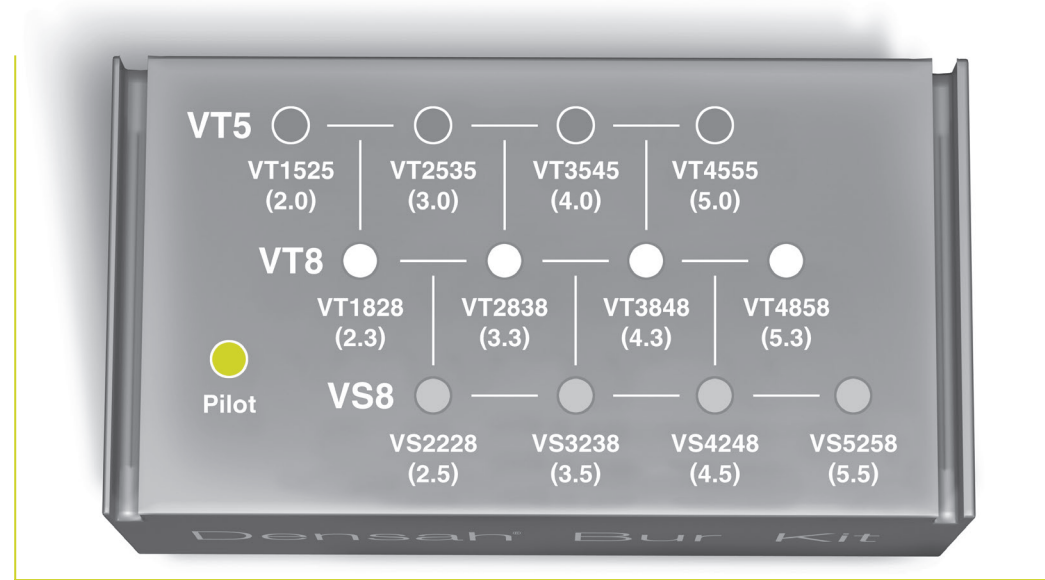

Use large block display to compare Bur identification system when using the schematic below for proper Bur usage

- VT5 Set
- VT8 Set
- VS8 Set

Densifying Mode CCW (800-1500) RPMs / Cutting Mode CW (800-1500) RPMs																	
Implant Direct®			SwishPlant™														
			Soft Bone							Hard Bone (Mandible)							
										In densifying mode make sure your osteotomy is 1.0 mm deeper than the actual implant final length. In extreme hard bone, utilize DAC (Densify After Cut) Protocol. Click here to view PDF.							
Geometry	Major Ø	Minor Ø	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Densah® Bur Block Display	Pilot	Bur 1	Bur 2	Bur 3	Bur 4	Bur 5	Bur 6	Densah® Bur Block Display
Taper	3.3	2.9	Pilot	VT1828 (2.3)	VS2228* (2.5)	—	—	—		Pilot	VT1828 (2.3)	VS2228* (2.5)	—	—	—	—	
Taper	4.1	3.5	Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—		Pilot	VT1828 (2.3)	VT2838 (3.3)	VS3238* (3.5)	—	—	—	
Taper	4.8	4.2	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545* (4.0)	—	—		Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848* (4.3)	—	
Taper	5.7	5.0	Pilot	VT1525 (2.0)	VT2535 (3.0)	VT3545 (4.0)	VT4555* (5.0)	—		Pilot	VT1525 (2.0)	VT2535 (3.0)	VT2838 (3.3)	VT3545 (4.0)	VT3848 (4.3)	VT4555* (5.0)	

*Denotes implant placement. *Clinician experience and judgement should be used in conjunction with the Densifying Reference Guide recommendation and suggested use protocols. W042 REV04